


BU-80

PERFORMANCE SIMPLICITY

DURABILITY

Auto-Clave Sterilizer BU-80 Series

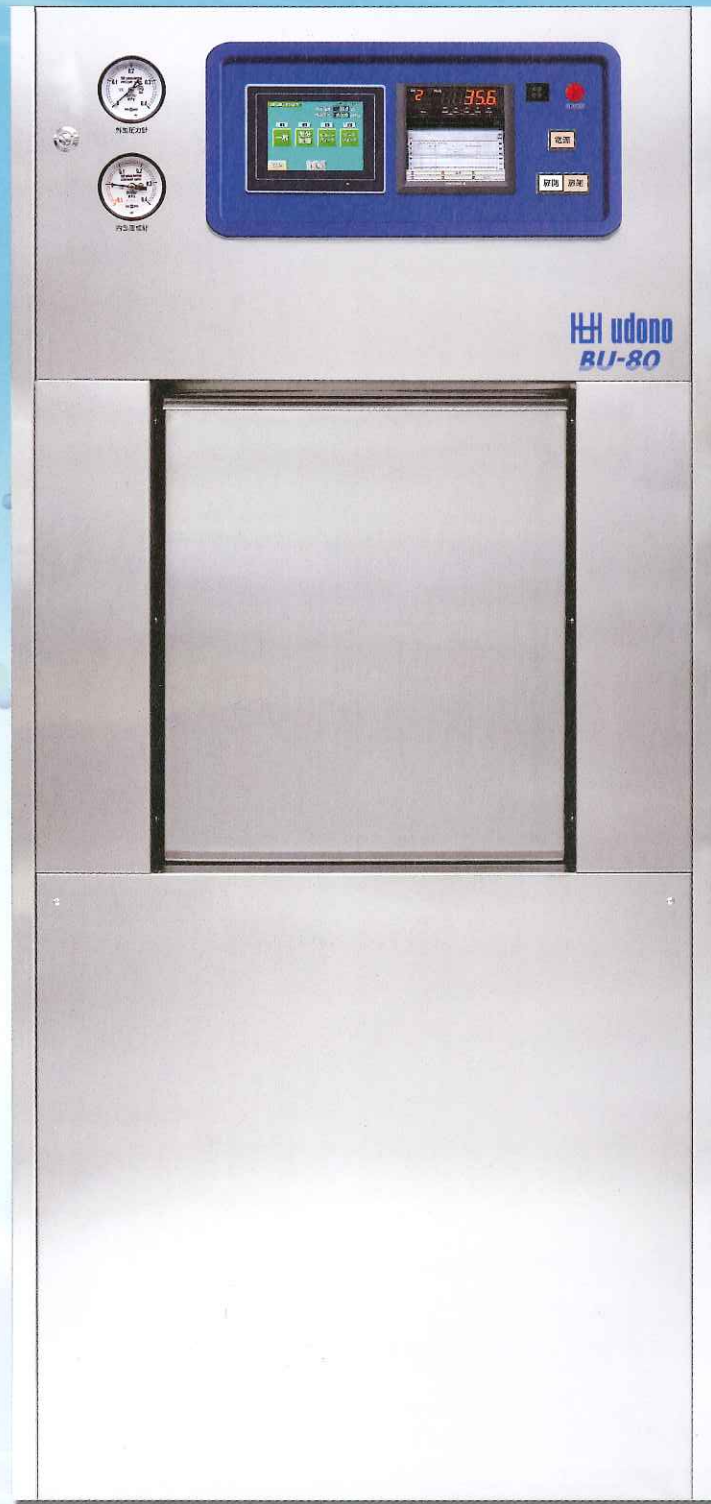


udono sterilization system



udono limited


**Udono brings this
high-performance Auto-Clave,
BU-80 series with
one-piece side wall construction, no weld
seams on the side-walls that could
increase wear completely covered corners
(for single door)
to enhance draining.**

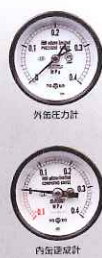


udono
BU-80




Friendly color touch panel screen

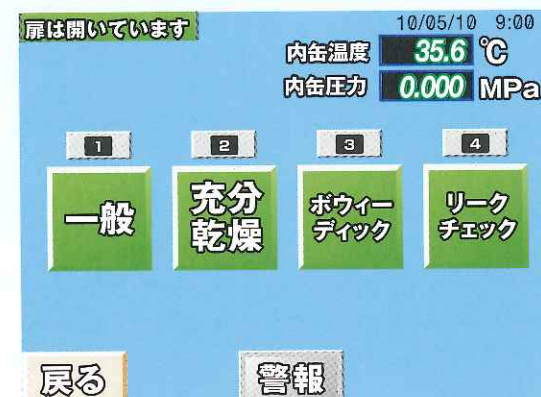


4 programs are available (General, Dry-well, BD and Leakage Test) which can be customized.

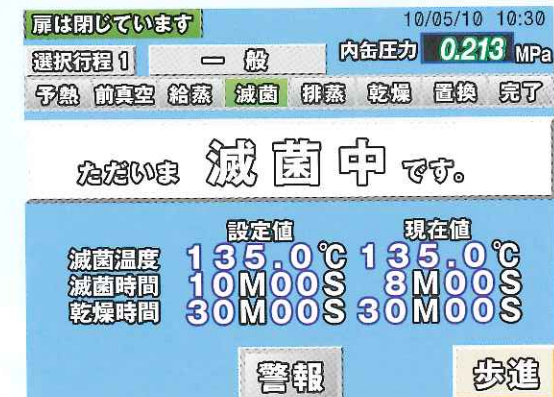

Intuitive color touch screen control system is welcomed. Process is shown on time.

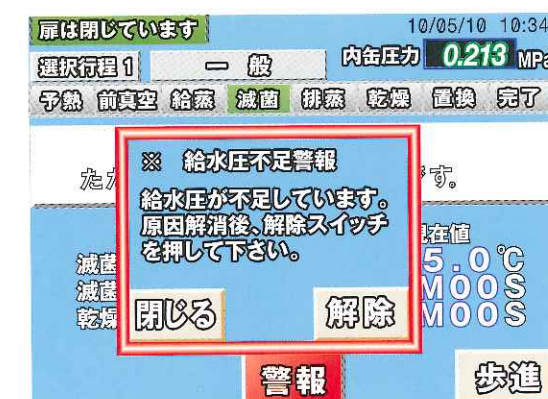


Trouble buzzer is equipped and what trouble happens is shown.

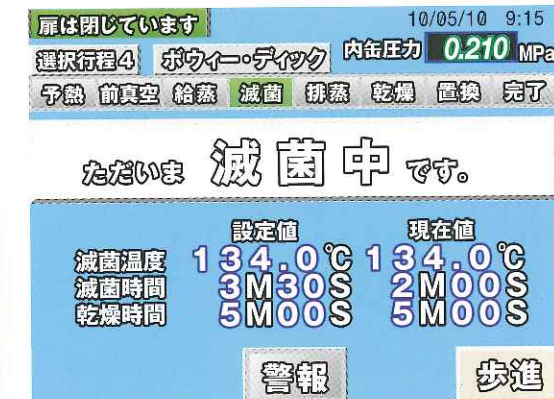

During BD, screen also shows the process on.
Screen can be showed in English

**BU-80 provides high-vacuum to
support "faster-dry" with
stronger chamber.**

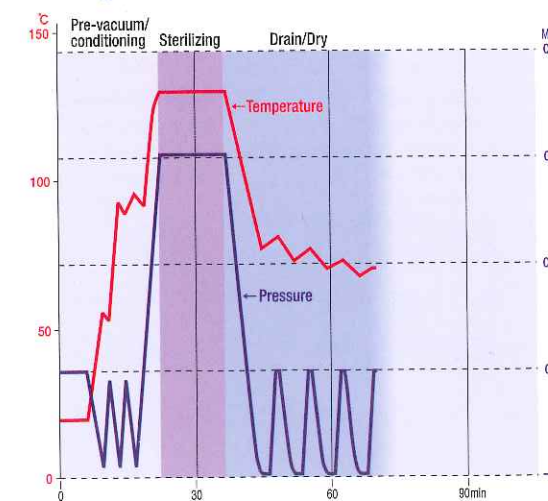


**Udono electric
heating clean
steam generator
is equipped as
standard.
(556, 559, 677, 6710
and 6713 only)**



Loading System (Option)



Shelves (T2 or T3)

Basket or Container is loaded in.
Shelves can be removed.

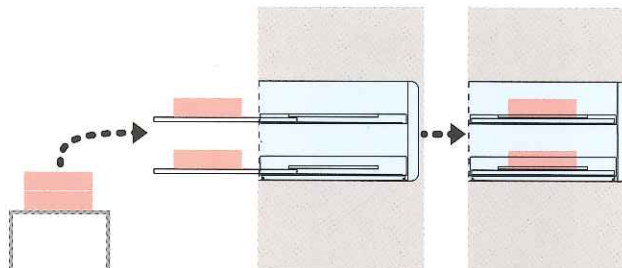


Trolley and loading rack (MT)

Basket is loaded on the rack and in together.
Trolley is used for carrying the rack with swivel casters.

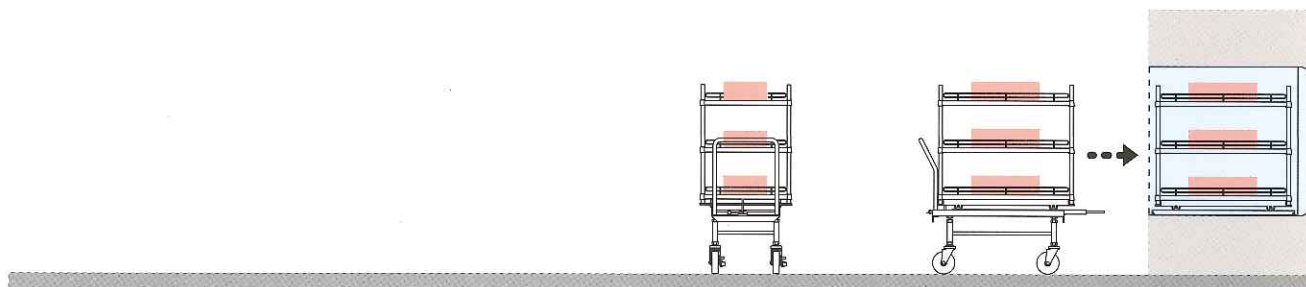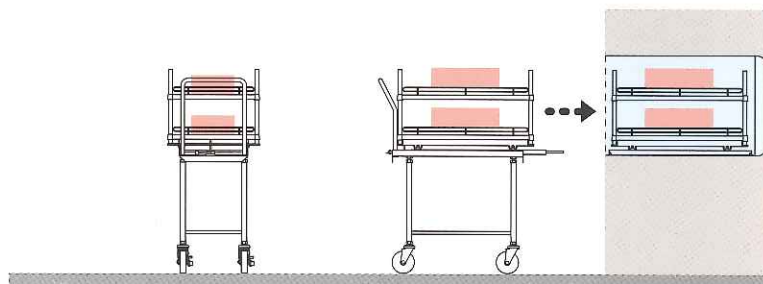


Folding rails and loading rack (SD)

Loading rack goes in and comes out.
Folding rails can be pulled down when loading rack comes out.

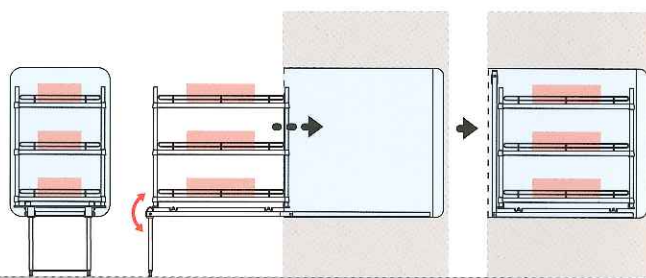


Semi-floor loading cart (AF)

The cart goes directly in.
This does not require recessed floor due of the wheels are pulled up.

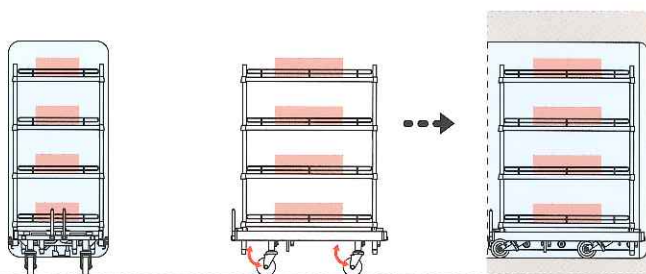

Floor loading cart (FL)

The cart goes in directly, however, recessed floor is required.

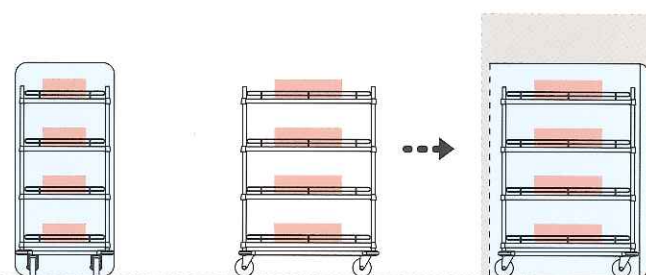

Door types

Slide door




BU-55 BU-67




BU-14

Horizontally or vertically can be selected. Vertical door matches to compact sized sterilizer. Horizontal door matches to larger sized sterilizer. The both door types can be operated by one touch.

Swing door



BU-R10



BU-R14

Door comes forward manually and is locked automatically.

Udono Clean Steam Generator

Udono clean steam generator provides “pure and clean steam” to protect instruments from getting rust but also it protects piping and chamber. High quality clean steam provided by Udono clean steam generator effectively protects “property of hospital = instruments” and “sterilizer” itself!

Specification
Slide door type

Specification Description	Chamber Size (mm)			Little	Door		Slide	Loadable basket or container	Overall Dimension (mm)				Steam Generator	Required Facilities					Power	
	Width	Height	Depth		Single	Double			Width	Height	Depth			Steam		Water Supply		Drain/ Exhaust Pipe dia.	AC200V 3 phases	Steam Generator AC200V 3 phases
											S	D		Single	Double	Pipe dia.	kg/hr			
BU-556-S	500	500	600	150	○	×	Vertically	4	850	1800	—	Yes	20A	15	6	15A	12kw 50A			
BU-559-*			900	230	○	○		6			1210			20						
BU-677-S	670	670	700	310	○	×		8	1100	—	25									
BU-6710-*			1000	440	○	○		12	1400	1310	35									
BU-6713-*			1300	580	○	○		16	1700	1610	45									
BU-6716-*			1600	710	○	○		20	2000	1910	55									
BU-107-S			670	1000	700	470		○	×	12	1075		—	35						
BU-1010-*					1000	670		○	○	18	1375		1410	55						
BU-1013-*	1300	870			○	○	24	1675	1710	70										
BU-1016-*	1600	1070			○	○	30	1975	2010	85										
BU-1410-*	1450	1000			970	○	○	24	1450	1410	75									
BU-1413-*		1300			1260	○	○	32	1750	1710	95									
BU-1416-*			1600	1550	○	○	40	2050	2010	115	12									

Swing door type

Specification Description	Chamber Size (mm)			Little	Door		Swing	Loadable basket or container	Overall Dimension (mm)			Steam Generator	Required Facilities					Power		
	Width	Height	Depth		Single	Double			Width	Height	Depth		Steam		Water Supply		Drain/ Exhaust Pipe dia.		AC200V 3 phases	
											S		D	Single	Double	Pipe dia.				kg/hr
BU-R107-S	670	1000	700	470	○	×	Right or Left	12	1070	1800	1000	—	NB	25A	35	8	40A	20A		
BU-R1010-*			1000	670	○	○		18			1300	1260			55					
BU-R1013-*			1300	870	○	○		24			1600	1560			70					
BU-R1016-*			1600	1070	○	○		30			1900	1860			85					
BU-R1410-*		1450	1000	970	○	○		24	1400	1290	75	32A		20A	10	50A	30A			
BU-R1413-*			1300	1260	○	○		32	1140	1850	1700								1590	95
BU-R1416-*			1600	1550	○	○		40		2000	1890								115	12

Description of Item #

BU-R1013-D-MT/ * *

Door Open
R: Swing
No mark: Slide

Item #

Door Type
S: Single
D: Passing through

Loading Type
T2,T3: Shelves
MT: Trolley and loading rack
SD: Folding rails and loading rack
AF: Semi-floor loading cart
FL: Floor loading cart

NB: No Steam Generator
NS: Normal Steam (556~6713)
No mark: Steam Generator (556~6713)



Head Office
2-1-9 Motoyokoyama-cho, Hachioji-shi, Tokyo 192-0063, Japan
Telephone +81-42-642-6301 E-mail sales@udono.com
Facsimile +81-42-644-8350 http://www.udono.com